

Separation of Lisinopril and its Organic Impurities per USP Monograph

Krishna Chaitanya Routhu¹, Swetha Kotikalapudi¹, and Bryan Tackett, PhD²

¹India Phenologix Lab, Phenomenex India, Hitech Defence and Aerospace Park Industrial Area, Mahadeva Kodigehalli, Holbi, Jala Taluka, Bengaluru 562149, India

²Phenomenex Inc., 411 Madrid Ave., Torrance, CA 90501 USA

Overview


Lisinopril is an orally active angiotensin-converting enzyme (ACE) inhibitor used for the treatment of hypertension, heart failure, and acute myocardial infarction. The development of a quick and efficient analysis of Lisinopril and its related organic impurities is of interest for generic drug manufacturers. In this application note, we report the separation of Lisinopril and its related organic impurities using a Luna™ 5 µm C8(2) column and a Kinetex™ 5 µm C8 column according to the USP monograph for Lisinopril.

System suitability per USP Monograph for the Lisinopril Assay is a symmetry factor no more than (NMT) 1.7 and a percent relative standard deviation (%RSD) of NMT 0.73 % for Lisinopril. System suitability per USP Monograph for the Lisinopril Organic Impurities is a %RSD of NMT 10.0 %, and a signal-to-noise (S/N) ratio no less than 10 for Lisinopril.

The results clearly show that the system suitability criteria (Symmetry, %RSD and S/N ratio) for both assay and organic impurities per the USP monograph for Lisinopril were met with both the Luna 5 µm C8(2) and Kinetex 5 µm C8 columns. While the data demonstrates that either column would be acceptable, there are certain advantages to be gained from the use of the core-shell Kinetex C8 column, specifically shorter run times and increased sensitivity (S/N), while the Luna C8(2) column provided slightly better peak shape.

All solutions were prepared as indicated in the USP Monograph for Lisinopril. USP Lisinopril RS (Catalog No. 1368609) was purchased from USP.

Figure 1. Lisinopril

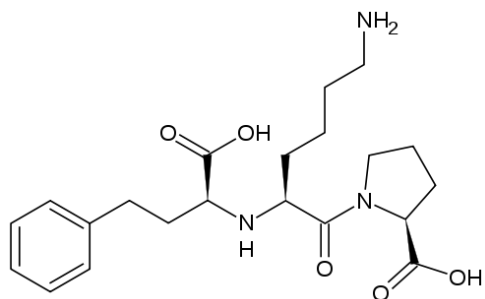

LC-UV Conditions - Assay

Column: Luna 5 µm C8(2) ([00G-4249-E0](#))

Kinetex 5 µm C8 ([00G-4608-E0](#))

Dimensions: 250 x 4.6 mm

Mobile Phase: Acetonitrile / Buffer (4:96, v/v)

Buffer: 2.76 g of Monobasic Sodium Phosphate was dissolved in 900 mL of water in a 1000 mL volumetric flask. Adjusted pH to 5.0 with 1N Sodium Hydroxide and diluted with water to volume.

Flow Rate: 1 mL/min (Isocratic)

Injection Volume: 20 µL

Temperature: 50 °C

Detector: UV @ 210 nm

System: Waters® ACQUITY Arc® HPLC

LC-UV Conditions - Organic Impurities

Column: Luna 5 µm C8(2) ([00G-4249-E0](#))

Kinetex 5 µm C8 ([00G-4608-E0](#))

Dimensions: 250 x 4.6 mm

Mobile Phase: A: Acetonitrile / Buffer (7:193, v/v)

B: Acetonitrile / Buffer (20:80, v/v)

Buffer: 3.53 g of Monobasic Sodium Phosphate Dihydrate is added to 1000 mL water and the pH is adjusted to 4.1 with Phosphoric Acid.

Gradient:	Time (min)	%B
	0	0
	35	40
	55	40
	60	0

Flow Rate: 1.8 mL/min

Injection Volume: 20 µL

Temperature: 45 °C

Detector: UV @ 210 nm

System: Waters ACQUITY Arc HPLC

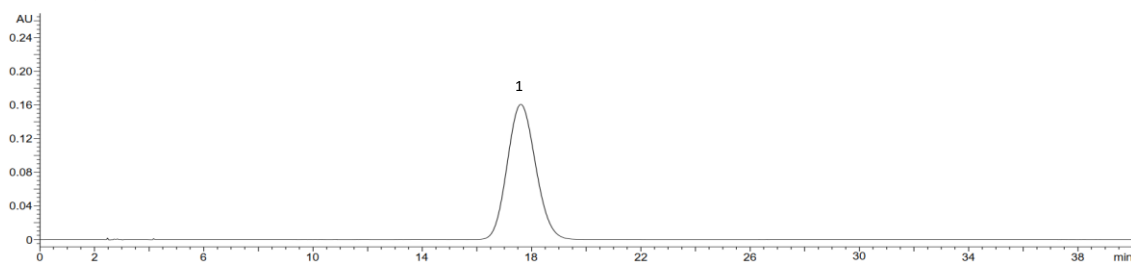


Table 1. Preparation of Solutions

Solution	Composition
Standard Solution – Assay	0.3 mg/mL of USP Lisinopril RS in water
Standard Solution – Organic Impurities	0.006 mg/mL of USP Lisinopril RS in Mobile Phase A
Sensitivity Solution – Organic Impurities	1.0 µg/mL of USP Lisinopril RS in Mobile Phase A from Standard Solution – Organic Impurities

Figure 2. Standard Solution – Assay

Luna™ 5 µm C8(2) Column




App ID 27240

Peak No.	Analyte	Retention Time (min)	Area	Area %RSD	Symmetry Factor
1	Lisinopril	17.64	11571239	0.04	1.11

N = 6 Injections

Kinetex™ 5 µm C8 Column

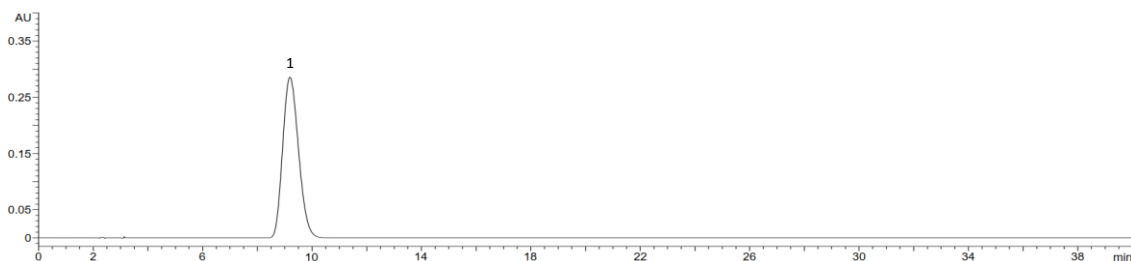


App ID 27241

Peak No.	Analyte	Retention Time (min)	Area	Area %RSD	Symmetry Factor
1	Lisinopril	9.20	11393656	0.35	1.19


N = 6 Injections



Figure 3. Standard Solution – Organic Impurities

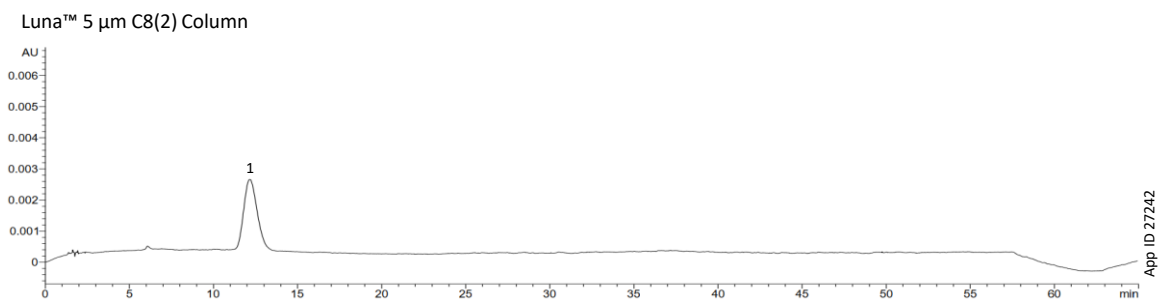


Peak No.	Analyte	Retention Time (min)	Area	Area %RSD	Symmetry Factor
1	Lisinopril	12.15	127254	0.98	1.23
N = 6 Injections					

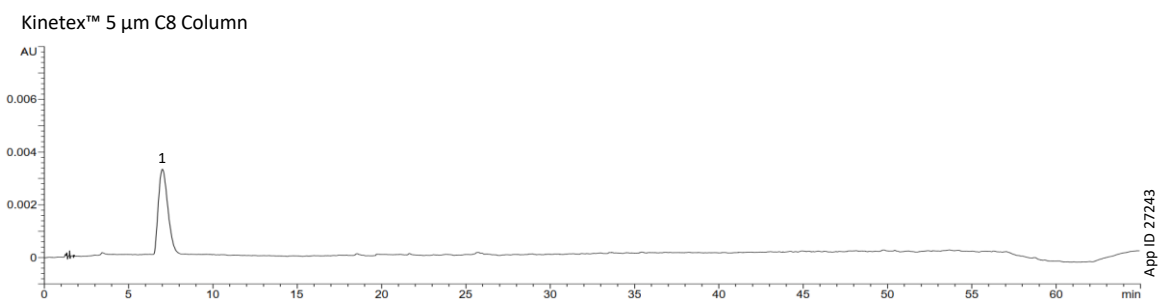

Peak No.	Analyte	Retention Time (min)	Area	Area %RSD	Symmetry Factor
1	Lisinopril	7.01	125862	0.73	1.38
N = 6 Injections					



Figure 4. Sensitivity Solution – Organic Impurities

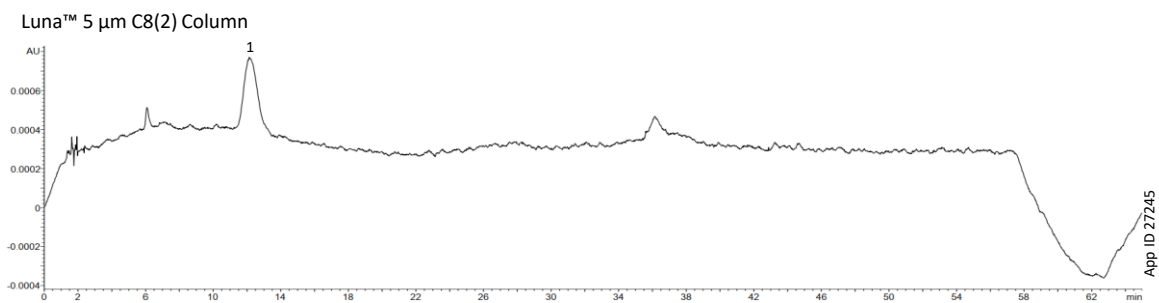

Peak No.	Analyte	Retention Time (min)	S/N Ratio
1	Lisinopril	12.14	35.15
N = 3 Injections			



Peak No.	Analyte	Retention Time (min)	S/N Ratio
1	Lisinopril	7.00	48.8
N = 3 Injections			



Need a different column size or sample preparation format?

No problem! We have a majority of our available dimensions up on www.phenomenex.com, but if you can't find what you need right away, our super helpful Technical Specialists can guide you to the solution via our online chat portal www.phenomenex.com/Chat.

Australia

t: +61 (0)2-9428-6444
auiinfo@phenomenex.com

Austria

t: +43 (0)1-319-1301
anfrage@phenomenex.com

Belgium

t: +32 (0)2 503 4015 (French)
t: +32 (0)2 511 8666 (Dutch)
beinfo@phenomenex.com

Canada

t: +1 (800) 543-3681
info@phenomenex.com

China

t: +86 400-606-8099
cninfo@phenomenex.com

Czech Republic

t: +420 272 017 077
cz-info@phenomenex.com

Denmark

t: +45 4824 8048
nordicinfo@phenomenex.com

Finland

t: +358 (0)9 4789 0063
nordicinfo@phenomenex.com

France

t: +33 (0)1 30 09 21 10
franceinfo@phenomenex.com

Germany

t: +49 (0)6021-58830-0
anfrage@phenomenex.com

Hong Kong

t: +852 6012 8162
hkinfo@phenomenex.com

India

t: +91 (0)40-3012 2400
indiainfo@phenomenex.com

Indonesia

t: +62 21 5019 9707
indoinfo@phenomenex.com

Ireland

t: +353 (0)1 247 5405
eireinfo@phenomenex.com

Italy

t: +39 051 6327511
italiainfo@phenomenex.com

Japan

t: +81 (0) 120-149-262
jpinfo@phenomenex.com

Luxembourg

t: +31 (0)30-2418700
nlinfo@phenomenex.com

Mexico

t: 01-800-844-5226
tecnicomx@phenomenex.com

The Netherlands

t: +31 (0)30-2418700
nlinfo@phenomenex.com

New Zealand

t: +64 (0)9-4780951
nzinfo@phenomenex.com

Norway

t: +47 810 02 005
nordicinfo@phenomenex.com

Poland

t: +48 22 104 21 72
pl-info@phenomenex.com

Portugal

t: +351 221 450 488
ptinfo@phenomenex.com

Singapore

t: +65 6559 4364
sginfo@phenomenex.com

Slovakia

t: +420 272 017 077
sk-info@phenomenex.com

Spain

t: +34 91-413-8613
espinfo@phenomenex.com

Sweden

t: +46 (0)8 611 6950
nordicinfo@phenomenex.com

Switzerland

t: +41 (0)61 692 20 20
swissinfo@phenomenex.com

Taiwan

t: +886 (0) 0801-49-1246
twinfo@phenomenex.com

Thailand

t: +66 (0) 2 566 0287
thaiinfo@phenomenex.com

United Kingdom

t: +44 (0)1625-501367
ukinfo@phenomenex.com

USA

t: +1 (310) 212-0555
www.phenomenex.com/chat

**🌐 All other countries/regions
Corporate Office USA**

t: +1 (310) 212-0555
www.phenomenex.com/chat

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country/region, contact Phenomenex USA, International Department at international@phenomenex.com

BE-HAPPY™
GUARANTEE

Your happiness is our mission. Take 45 days to try our products. If you are not happy, we'll make it right.

www.phenomenex.com/behappy

Terms and Conditions

Subject to Phenomenex Standard Terms and Conditions, which may be viewed at www.phenomenex.com/phx-terms-and-conditions-of-sales.

Trademarks

Luna, Kinetex, and BE-HAPPY are trademarks of Phenomenex. Waters and ACQUITY Arc are registered trademarks of Waters Technologies Corporation.

Disclaimer

Comparative separations may not be representative of all applications. Phenomenex is in no way affiliated with Waters Technologies Corporation.

FOR RESEARCH USE ONLY. Not for use in clinical diagnostic procedures.

© 2022 Phenomenex, Inc. All rights reserved.

