

USP Albuterol Inhalation Solution: Assay on Kinetex® 2.6 µm C18 and Organic Impurities on Kinetex® 1.7 µm XB-C8

Brittany LeBlanc and Phil Koerner
Phenomenex, Inc., 411 Madrid Ave., Torrance, CA 90501 USA

Overview

Albuterol Inhalation Solution is an isotonic sterile solution of albuterol sulfate; albuterol (aka Salbutamol) is a bronchodilator used to treat or prevent bronchospasm in people with reversible obstructive airway disease. This application note for the LC-UV assay and organic impurities for Albuterol inhalation solution is based on the proposed USP monograph published in PF45(3). The mobile phase and standard solution for the assay was prepared in accordance with the proposed USP monograph, which is not yet official. The proposed monograph was validated using Kinetex 2.6 µm C18 for the assay, and validated using Kinetex 1.7 µm XB-C18 for the organic impurities.

USP reference standards were used to prepare the standard solutions and system suitability solution, and mobile phase was prepared per the USP draft monograph for Albuterol inhalation solution.


LC-UV Conditions (Assay)

Column: Kinetex 2.6 µm C18 ([00F-4462-E0](#))
Dimension: 150 x 4.6 mm
Mobile Phase: A = 3.4 g/L of monobasic potassium phosphate and 1.1 g of sodium 1-heptanesulfonate in water. Adjust with phosphoric acid to pH 2.1.
 B = Acetonitrile

Gradient: Time (min)	% B
0	15
6	40
7.5	40
7.6	15
11	15

Flow Rate: 0.75 mL/min
Temperature: 37 °C
Detector: UV @ 210 nm
Injection: 10 µL
System: Agilent® 1290 Infinity II (binary)
Sample: As indicated in chromatogram

Standard Solution for Assay on Kinetex 2.6 µm C18

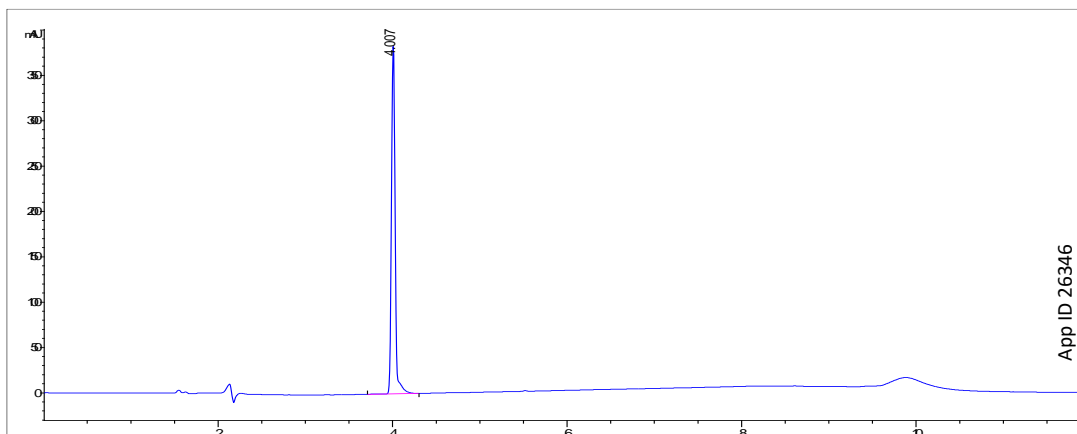

Standard Solution
 (0.1 mg/mL of USP Albuterol Sulfate RS (equivalent to 0.08 mg/mL of albuterol) in 0.01N HCl)

Table 1. Summary of Results, Assay (Standard Solution)

Analyte	Retention Time, min	Symmetry Factor (NMT 1.7)	% RSD (n = 6) (NMT 1.0 %)
Albuterol sulfate	4.007	1.08	0.16




LC-UV Conditions (Organic Impurities)

Column: Kinetex 1.7 µm XB-C18 ([00F-4498-AN](#))
Dimension: 150 x 2.1 mm
Mobile Phase: A = 3.4 g/L of monobasic potassium phosphate and 1.1 g of sodium 1-heptanesulfonate in water. Adjust with phosphoric acid to pH 2.1.
 B = Acetonitrile

Gradient: Time (min)	% B
0	5
2.5	7.5
5	15
18	19.5
26	36
26.5	50
27.5	50
27.6	5
34	5

Flow Rate: 0.35 mL/min
Temperature: 37 °C
Detector: UV @ 210 nm
Injection: 10 µL
System: Agilent® 1290 Infinity II (binary)
Sample: As indicated in chromatograms


Standard Solution for Organic Impurities on Kinetex 1.7 µm XB-C18

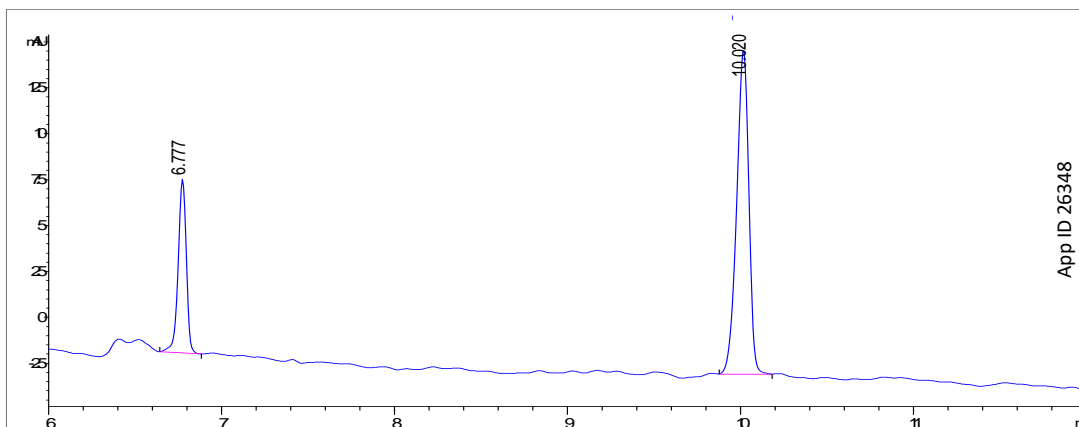

Standard Solution
 (1.25 µg/mL each of USP Albuterol Sulfate RS and USP Levalbuterol Related Compound D RS in 0.01N HCl)

App ID 26348

System Suitability Solution for Organic Impurities on Kinetex 1.7 µm XB-C18

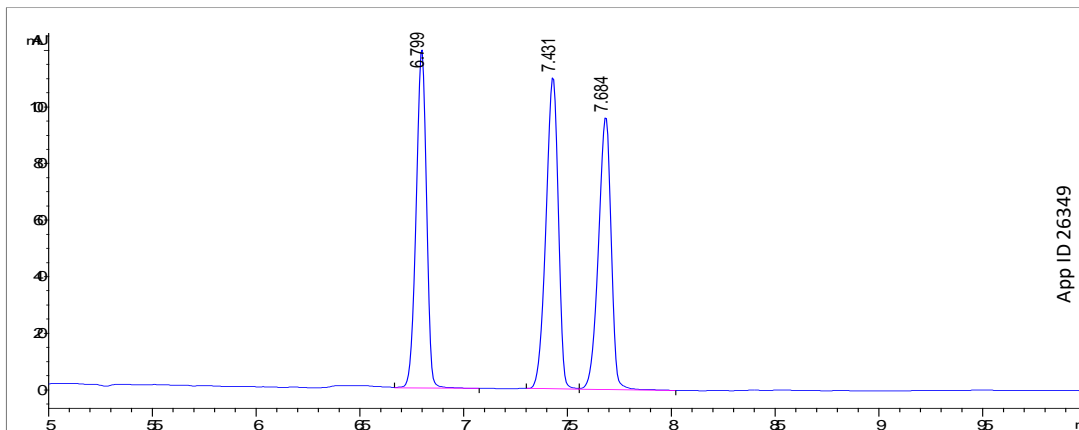

System Suitability Solution
 (0.05 mg/mL each of USP Albuterol Sulfate RS, USP Albuterol Related Compound I RS, and USP Levalbuterol Related Compound H RS in 0.01N HCl)

App ID 26349



Table 2. Summary of Results, Organic Impurities (Standard Solution)

Analyte	Retention Time, min	Symmetry Factor (NMT 2.0)	% RSD (n = 6) (NMT 5.0 %)	Signal-to-Noise (NLT 10)
Albuterol sulfate	6.777	1.03	1.3	15
Levalbuterol Related Compound D	10.020	0.96	1.7	50

Table 3. Summary of Results, Organic Impurities (System Suitability Solution)

Analyte	Retention Time, min	Resolution (NLT 2.0)
Albuterol sulfate	6.805	
Albuterol Related Compound I	7.434	5.9
Albuterol Related Compound H	7.687	2.2

Order Information

Kinetex 2.6 μm Analytical Columns (mm)

Phases 150 x 4.6

C18 [00F-4462-E0](#)

Kinetex 1.7 μm Analytical Columns (mm)

Phases 150 x 2.1

XB-C18 [00F-4498-AN](#)



Need a different column size or sample preparation format?

No problem! We have a majority of our available dimensions up on www.phenomenex.com, but if you can't find what you need right away, our super helpful Technical Specialists can guide you to the solution via our online chat portal www.phenomenex.com/LiveChat.

Australia

t: +61 (0)2-9428-6444
auinfo@phenomenex.com

Austria

t: +43 (0)1-319-1301
anfrage@phenomenex.com

Belgium

t: +32 (0)2 503 4015 (French)
t: +32 (0)2 511 8666 (Dutch)
beinfo@phenomenex.com

Canada

t: +1 (800) 543-3681
info@phenomenex.com

China

t: +86 400-606-8099
cninfo@phenomenex.com

Czech Republic

t: +420 272 017 077
cz-info@phenomenex.com

Denmark

t: +45 4824 8048
nordicinfo@phenomenex.com

Finland

t: +358 (0)9 4789 0063
nordicinfo@phenomenex.com

France

t: +33 (0)1 30 09 21 10
franceinfo@phenomenex.com

Germany

t: +49 (0)6021-58830-0
anfrage@phenomenex.com

Hong Kong

t: +852 6012 8162
hkinfo@phenomenex.com

India

t: +91 (0)40-3012 2400
indiainfo@phenomenex.com

Ireland

t: +353 (0)1 247 5405
eirinfo@phenomenex.com

Italy

t: +39 051 6327511
italiainfo@phenomenex.com

Luxembourg

t: +31 (0)30-2418700
nlinfo@phenomenex.com

Mexico

t: 01-800-844-5226
tecnicomx@phenomenex.com

The Netherlands

t: +31 (0)30-2418700
nlinfo@phenomenex.com

New Zealand

t: +64 (0)9-4780951
nzinfo@phenomenex.com

Norway

t: +47 810 02 005
nordicinfo@phenomenex.com

Poland

t: +48 22 104 21 72
pl-info@phenomenex.com

Portugal

t: +351 221 450 488
ptinfo@phenomenex.com

Singapore

t: +65 800-852-3944
sginfo@phenomenex.com

Slovakia

t: +420 272 017 077
sk-info@phenomenex.com

Spain

t: +34 91-413-8613
espinfo@phenomenex.com

Sweden

t: +46 (0)8 611 6950
nordicinfo@phenomenex.com

Switzerland

t: +41 (0)61 692 20 20
swissinfo@phenomenex.com

Taiwan

t: +886 (0) 0801-49-1246
twinfo@phenomenex.com

Thailand

t: +66 (0) 2 566 0287
thaiinfo@phenomenex.com

United Kingdom

t: +44 (0)1625-501367
ukinfo@phenomenex.com

USA

t: +1 (310) 212-0555
info@phenomenex.com

Ⓢ All other countries/regions

Corporate Office USA
t: +1 (310) 212-0555
info@phenomenex.com

www.phenomenex.com

Phenomenex products are available worldwide. For the distributor in your country/region, contact Phenomenex USA, International Department at international@phenomenex.com

BE-HAPPY™
GUARANTEE

Your happiness is our mission. Take 45 days to try our products. If you are not happy, we'll make it right.

www.phenomenex.com/behappy

Terms and Conditions

Subject to Phenomenex Standard Terms and Conditions, which may be viewed at www.phenomenex.com/TermsAndConditions.

Trademarks

Kinetex is a registered trademark and BE-HAPPY is a trademark of Phenomenex. Agilent is a registered trademark of Agilent Technologies.

Disclaimer

Comparative separations may not be representative of all applications.

Phenomenex is in no way affiliated with Agilent Technologies.

FOR RESEARCH USE ONLY. Not for use in clinical diagnostic procedures.

© 2021 Phenomenex, Inc. All rights reserved.

